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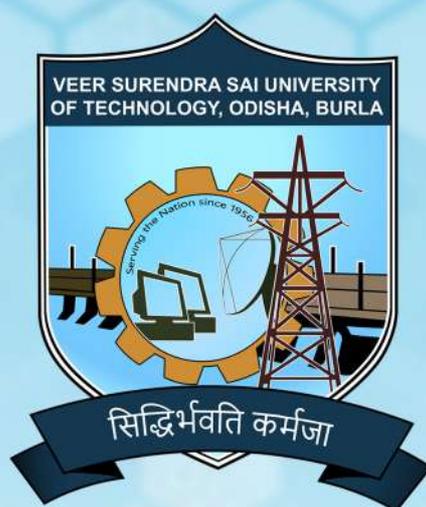
Regd. No..... Sem.....

DEPARTMENT OF CHEMISTRY

VEER SURENDRA SAI UNIVERSITY OF TECHNOLOGY

Siddhi Vihar, Burla-768018, Sambalpur, Odisha, India


www.vssut.ac.in

LABORATORY

OBSERVATION BOOK

Name

Regd. No

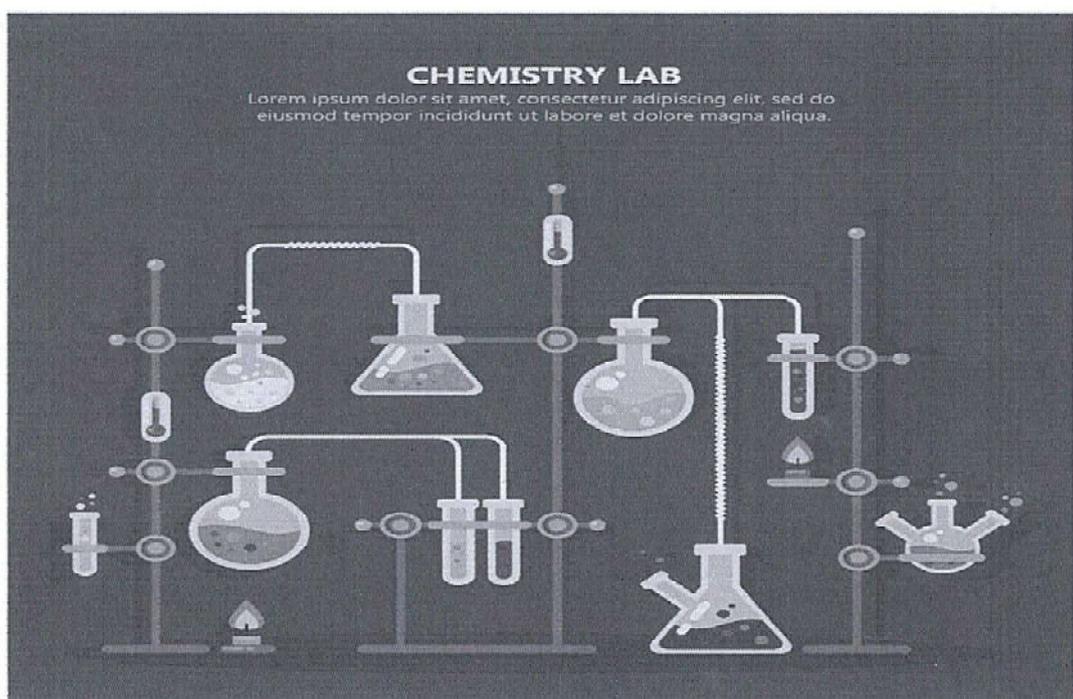DEPARTMENT OF CHEMISTRY

VEER SURENDRA SAI UNIVERSITY OF TECHNOLOGY, BURLA
ODISHA - 768018

LABORATORY PRACTICAL COPY
LABORATORY OBSERVATION BOOK

Name of the Student :.....

Roll No. :

Semester :

Academic Session :.....

DEPARTMENT OF CHEMISTRY
VEER SURENDRA SAI UNIVERSITY OF TECHNOLOGY, BURLA
SAMBALPUR, ODISHA - 768018

TABLE OF CONTENTS

Sl. No.	Title of Experiment	Page No.
1	Experiment - 1	1
2	Experiment - 2	3
3	Experiment - 3	6
4	Experiment - 4	8
5	Experiment - 5	9
6	Experiment - 6	11
7	Experiment - 7	13
8	Experiment - 8	16
9	Experiment - 9	19
10	Experiment - 10	22

CHEMISTRY LABORATORY
ODD/EVEN SEMESTER 20.....
LOGBOOK MARK SHEET

Name :

Regd. No. :

Branch :

Section :

Semester :

TABLE OF CONTENTS			
Sl. No.	Title of Experiment	Marks obtained	Signature
1	 / 10	
2	 / 10	
3	 / 10	
4	 / 10	
5	 / 10	
6	 / 10	
7	 / 10	
8	 / 10	
9	 / 10	
10	 / 10	
	Total Marks obtained / 100	

Unauthorized absences will be credited with a mark of 0/10. Your completed laboratory logbook (with this sheet firmly stuck in and fully completed) must be submitted to your instructor at the end of all experiments. Your instructor will make an assessment and transmit the mark

Experiment :1

Aim of the Experiment : To determine the total amount of Sodium Carbonate and Sodium Hydroxide present in the supplied sample of Water.

Apparatus Required : Burette, Pipette, Conical Flask, Ice Bath.



Chemicals Required : Standard Acid, Methyl Orange, Phenolphthalein.

Theory : When a mixture of a base and a carbonate is titrated with a standard acid, then the neutralization occurs in two steps.

(i) In the first step, NaOH is completely neutralized and Na_2CO_3 is partially neutralized.

(ii) In the second step, neutralization of Na_2CO_3 is completed.

Chemical Reaction :

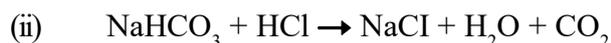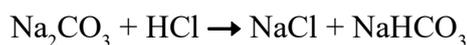Procedure:

Pipette out 10 ml of supplied water sample into a conical flask. Add two drops of phenolphthalein indicator and cool it in an ice bath till it gets sufficiently cooled. Then titrate against standard acid (HCl) till the pink coloured solution turned into a colourless solution. Note the burette reading. The difference between the final and initial burette reading gives the volume V_1 . To this titrated solution, add 2 drops of methyl orange indicator and continue the titration till the colour changes from yellow to red. Again note down the final reading and the required volume of standard acid, V_2 . Repeat the above process three times to get the concordant readings.

Tabulation :

Number of Observation	Volume of Water Sample in ml	Burette Reading in ml			Difference in ml		Remarks
		I.B.R. (a)	F.B.R. (1 st) (b)	F.B.R. (2 nd) (c)	(b-a) V_1	(c-a) V_2	
1							Rough
2							concordant
3							
4							

I.B.R. Initial Burette Reading



F.B.R. Final Burette Reading

Calculations:




Subtracting equation (I) from equation (II)

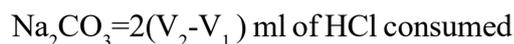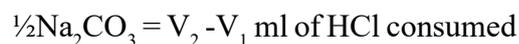Similarly, the volume of Acid required to neutralize whole of NaOH = $V_2 - 2(V_2 - V_1)$ ml

Amount of Na_2CO_3 :

We know that

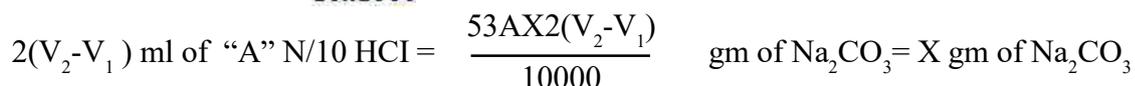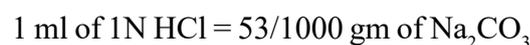
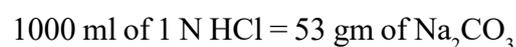
10 ml of supplied water sample contains X gm of Na_2CO_3 ,

1000 ml of Supplied water sample contains 100X gm/L of Na_2CO_3 ,

Amount of NaOH:

We know that

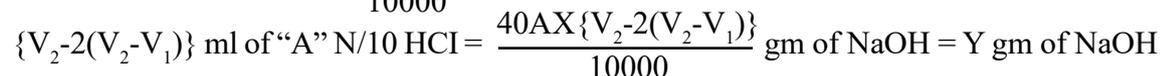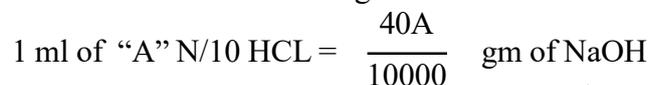10 ml of supplied water sample contains Y gm of NaOH

1000 ml of Supplied water sample contains 100Y gm/L of NaOH.

Conclusion : The supplied water sample contains 100X gm/L of sodium carbonate and 100Y gm/L of Sodium hydroxide.

Experiment : 2

Aim of the Experiment : To determine the hardness of water by Ethylene Diamine Tetraacetic Acid (EDTA) method.

Apparatus Required : Burette, Pipette, Conical Flask, Measuring Flask, Wash Bottle

Chemicals Required:

(i) Standard disodium-EDTA Solution (ii) Eriochrome black-T indicator (EBT) [Sodium-1-(1-hydroxy-2-naphthyl azo)-6-nitro-2-naphthol-4-sulphonate (II)] (iii) Buffer solution ($\text{NH}_4\text{Cl} + \text{NH}_4\text{OH}$)

Theory : Hardness is the property of water, which prevents the lathering of soap. Originally, it was defined as soap consuming capacity of the water sample. The hardness of water is generally due to the presence of certain salts of Ca^{2+} , Mg^{2+} and other heavy metal ions like Al^{3+} , Fe^{3+} and Mn^{2+} dissolved in it. A sample of hard water, when treated with soap (sodium or potassium salt of higher fatty acid like oleic acid, palmitic acid or stearic acid) does not produce foam but on the other hand forms insoluble white scum or precipitate, which do not possess any detergent action. This is due to the formation of insoluble soaps of calcium and magnesium.

Type of Hardness:

(a) Temporary Hardness : It is caused by the presence of dissolved bicarbonates of calcium, magnesium and other heavy metals and the carbonates of iron. Temporary hardness can be removed by boiling of water as bicarbonate salts are decomposed producing insoluble carbonates or hydroxides which are deposited as a crust at the bottom of the vessel.

(b) Permanent Hardness : It is due to the presence of dissolved chlorides and sulphates of calcium, magnesium, iron and other heavy metals. Hence, the salts responsible for permanent hardness are CaCl_2 , MgCl_2 , CaSO_4 , MgSO_4 , FeSO_4 , $\text{Al}_2(\text{SO}_4)_3$, etc. Hardness of water is generally expressed in terms of equivalent amount of CaCO_3 . The reason for choosing CaCO_3 as the standard for reporting hardness of water is actually its molecular weight which is 100. This makes the calculation easier. Moreover, CaCO_3 is the most insoluble salt that can be precipitated in the water treatment.

Units of Hardness:

(i) Parts per million (PPM) : It is defined as the number parts by weight of calcium carbonate present per million (10^6) parts by weight of water i.e. 1 PPM = 1 part of CaCO_3 equivalent hardness in 10^6 parts of water.




(ii) Milligrams per Liter (mg/L) : It is defined as the number of milligrams of calcium carbonate present in one litre of water.

It should be remembered that $1\text{mg/L} = 1\text{PPM}$ for water.

Principle : Disodium salt of EDTA (a chelating ligand) when treated with hard water reacts with it and consumes all the hardness producing cations in alkaline medium (pH 8-10). The excess of reagent is indicated by the reaction of the indicator Eriochrome black-T, which colours the solution from wine red to clear blue that is the end point of this complexometric titration.

Chemical Reactions :

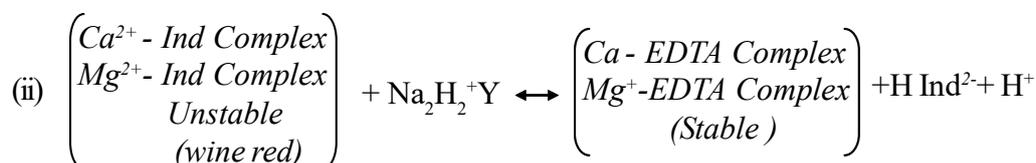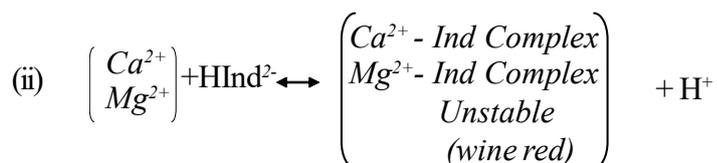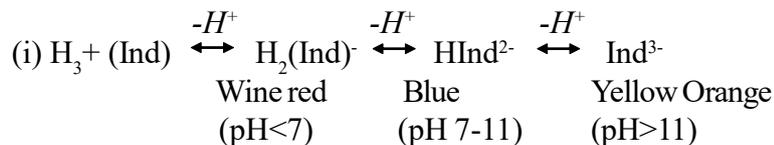Procedure:

(i) Determination of Total Hardness:

Pipette out 40 ml of hard water sample into a 100 ml of conical flask. Add 2 ml of buffer solution and 2-3 drops of Ericrome black-T indicator. Then titrate the solution with M/100 EDTA solution from the burette until colour changes from wine red to blue colour at the end point. Note the titrate value till you get the concordant reading. Let the volume is V_2 ml.

Tabulation: 01

No of Observations	Amount of supplied water sample in ml	Initial Burette reading in ml	Final Burette reading in ml	Difference in ml	Remarks
1					Rough
2					Concordant
3					
4					

(ii) Determination of Temporary and Permanent Hardness : Take 250 ml of hard water sample in a large beaker and boil gently for about 1 hour (till all the bicarbonates are decomposed into insoluble $\text{CaCO}_3 + \text{Mg}(\text{OH})_2$). Cool, filter into a 250 ml measuring flask and make the volume up to the mark. Take 50 ml of this solution and proceed in the same way as above. The volume of EDTA used corresponds to the permanent hardness in the water sample. Let the volume is V_3 ml.

Temporary hardness is calculated by subtracting permanent hardness from total hardness.

**Tabulations: 02**

No of Observations	Amount of supplied water sample in ml	Initial Burette reading in ml	Final Burette reading in ml	Difference in ml	Remarks
1					Rough
2					Concordant
3					
4					

Calculation :**Total Hardness:**

1000 ml of 1 M EDTA = 100 gm of CaCO_3

1 ml of 1 M EDTA = 100 mg of CaCO_3

1 ml of "a" M/100 EDTA = "a" mg of CaCO_3

V_2 ml of "a" M/100 EDTA = "a" V_2 mg of CaCO_3 = "A" mg of CaCO_3

Now 40 ml water sample contains = "A" mg of CaCO_3

1000 ml or 1 liter of water sample contains = $\frac{A}{40} \times 1000 = 25A$ mg/L = 25A ppm total

hardness in terms of CaCO_3 equivalent hardness.

Permanent Hardness:

1 ml of 1 "a" M/100 EDTA = "a" mg of CaCO_3

V_3 ml of "a" M/100 EDTA = "a" V_3 mg of CaCO_3 = "B" mg of CaCO_3

Now 40 ml boiled water sample contains = "B" mg of CaCO_3

1000 ml or 1 Liter of boiled water sample contains = $\frac{B}{40} \times 1000 = 25B$ mg/L = 25B ppm

permanent hardness in terms of CaCO_3 equivalent hardness.

So temporary hardness is equal to (25A-25B) ppm in terms of CaCO_3 equivalent hardness.

Conclusion : The supplied water sample contains 25A ppm of total hardness, 25B ppm permanent hardness and (25A-25B) ppm of temporary hardness in terms of CaCO_3 equivalent hardness.

Experiment : 3


Aim of the Experiment : Standardization of Potassium Permanganate (KMnO_4) using Oxalic Acid

Apparatus Required : Burette, Pipette, Conical Flask

Chemicals Required : KMnO_4 Solution, Oxalic Acid (N/10), Dilute H_2SO_4

Theory : Potassium permanganate is a secondary standard solution since it contains traces of manganese dioxide as impurity. Moreover, potassium permanganate is a strong oxidizing agent and prone to decompose in presence of water, air and light. As a result, the strength of the potassium permanganate does not remain fixed for a longer time or period. Therefore, the strength of the permanganate solution is determined prior to its use. Potassium permanganate acts as a strong oxidizing agent in acidic medium and therefore it can be titrated against a reducing agent like standard oxalic acid at 60-70 °C in presence of sulphuric acid. Oxalic acid can be decomposed on high-temperature heating before titration. During the reaction, MnO_4^- is reduced to Mn^{2+} and oxalate ions is oxidized to CO_2 . In this titration, potassium permanganate acts as a self-indicator because it has the capability to change colour with change in oxidation state (Mn^{7+} is pink in colour, whereas Mn^{2+} is colorless in water medium). During titration potassium permanganate is reduced by oxalic acid to colourless manganous ions and when the whole reducing agent is consumed, the addition of even one extra drop of KMnO_4 imparts light pink colouration to the solution, which is the end point of titration.

Chemical Reaction:

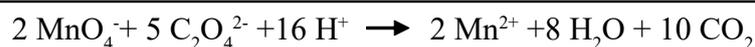Procedure:

Pipette out 10 ml of given N/10 oxalic acid into a 250 ml of conical flask. Add half test tubes of dilute sulphuric acid and heat up to 50-60 °C. Then, titrate the solution against KMnO_4 solution until to get a permanent light pink color. Repeat the titration till concordant readings are obtained.

Tabulation

No of Observations	Volume of N/10 Oxalic acid	I.B.R. (ml)	F.B.R (ml)	Difference in ml	Remarks
1					Rough
2					Concordant
3					
4					

Calculation :

We know; $S_1V_1 = S_2V_2$

Where, S_1 - strength of Oxalic acid (N/10)

V_1 - Volume of Oxalic acid (10 ml)

S_2 - Strength of $KMnO_4$

V_2 - Volume of $KMnO_4$ consumed

$$S_2 = \frac{S_1 \times V_1}{V_2}$$

Conclusion : Strength of supplied potassium permanganate isN.

Experiment : 4


Aim of the Experiment : Determination of Ferrous ion in Mohr's salt $[\text{FeSO}_4(\text{NH}_4)_2\text{SO}_4 \cdot 6\text{H}_2\text{O}]$ by using standardized potassium permanganate.

Apparatus Required : Burette, Pipette, Conical Flask, Beaker, Measuring Cylinder, Wash Bottle, Test Tube.

Chemicals Required : (i) Mohr's salt solution, (ii) KMnO_4 solution, (iii) Dilute H_2SO_4

Theory: The amount of ferrous ion present in Mohr's salt solution is determined by titrating against potassium permanganate solution in acidic medium. Here, potassium permanganate acts as a self-indicator. The end point is determined by the appearance of permanent light pink color. During the reaction Fe^{2+} is oxidized to Fe^{3+} and MnO_4^- is reduced to Mn^{2+}

Chemical Reactions:

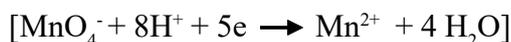Procedure : Pipette out 10 ml of supplied Mohr's salt solution into a 250 ml of conical flask. Add half test tubes of dilute sulphuric acid and titrate the solution against KMnO_4 solution until to get permanent light pink color. Repeat the titration till concordant readings are obtained.

No of Observations	Volume of N/10 Oxalic acid	I.B.R. (ml)	F.B.R. (ml)	Difference in ml	Remarks
1					Rough
2					Concordant
3					
4					

Calculation : Let the volume of KMnO_4 solution consumed by Mohr's salt solution is = V ml

We know 1000 ml of 1 N $\text{KMnO}_4 = 55.84$ gm of Fe

$$1 \text{ ml of } 1 \text{ N } \text{KMnO}_4 = \frac{55.85}{1000} \text{ gm of Fe}$$

$$V \text{ ml of } \frac{x}{10} \text{ N } \text{KMnO}_4 = \frac{55.85}{1000} \times V \times \frac{x}{10} \text{ gm of Fe} = 'A' \text{ g of Fe}$$

So 10 ml of Mohr's salt solution contains "A" gm of Fe

1000 ml of Supplied Mohr's salt solution contains 100 A gm/L of Fe.

Conclusion : Supplied Mohr's salt solution contains..... gm/L of Fe.

Experiment : 5

Aim of the Experiment : To determine the amount of Dissolved Oxygen present in the supplied sample of water.

Apparatus Required : Burette, Stopper Bottle (500 ml), Measuring Cylinder (250 ml), Test Tube.

Chemicals required :


(i) Manganous sulphate solution (48%), (ii) Alkaline potassium iodide solution, (iii) Concentrated sulphuric acid, (iv) Standard sodium thiosulphate solution, and (v) Freshly prepared starch solution.

Theory:

Dissolved oxygen is an important water quality parameter. It is the amount of oxygen required for complete oxidation of organic matter present in water. Water usually contains about 8 ml of dissolved oxygen per liter of water at room temperature. Determining dissolved oxygen is essential in boiler feed water because it can cause corrosion of the boiler. Dissolved oxygen in water is determined by Winkler's method.

Determination of dissolved oxygen in water is based on the oxidation of potassium iodide by dissolved oxygen. The iodometric method (Winkler) is based on the fact that in alkaline medium, dissolved oxygen oxidizes Mn^{+2} to Mn^{+4} , which in acidic medium oxidizes iodide to free iodine. The amount of iodine released is titrated against standard sodium thiosulphate solution.



Since dissolved oxygen in water is in the molecular state and is not capable of reacting with potassium iodide solution, oxygen carrier (such as Manganous hydroxide) is used to bring about the reaction between KI and oxygen.

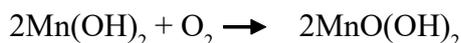(Oxygen carrier) (Basic Manganic Oxide)

(Nascent Oxygen)

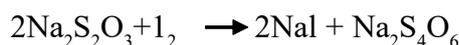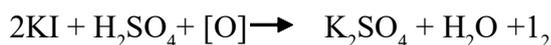Starch + I_2 Starch Iodide complex (Blue colour)

Procedure:

Take 250 ml of water sample in a stoppered bottle leaving behind around 10 ml of empty space at the top. (Avoiding as far as possible to contact with air). Immediately add 5 ml of manganous sulphate solution

followed by 5 ml of alkaline potassium iodide solution with the help of a pipette so that the solution is delivered well below the surface of water. (The alkaline KI solution contains sodium azide, which takes care of any oxidizing agent present, that are likely to liberate iodine from KI). The bottle is shaken well for several times to ensure proper mixing and the precipitate formed is allowed to settle down for a few minutes. Then 3 ml of concentrated H_2SO_4 is added dropwise to the solution with shaking till the precipitate is getting dissolved and giving a clear yellow color solution. Keep the solution in dark for about 5 minutes inside the cupboard. Then, pipette out 20 ml of the solution and titrate the solution against standard sodium thiosulphate solution from the burette. When the color of the solution is very light yellow (straw color), add 1 ml of freshly prepared starch indicator and continue the titration till the blue color just disappears.

No of Observations	Volume of N/10 Oxalic acid	Burette Reading (ml)		Difference in ml	Remarks
		I.B.R. (ml)	F.B.R (ml)		
1					Rough
2					Concordant
3					
4					

Calculations:

1000 ml of $1N Na_2S_2O_3$ solution equiv 8 gm of Oxygen

1 ml of $1N Na_2SO_3$ solution equiv 8 mg of Oxygen

V ml of "a" N / 10 $Na_2S_2O_3$ solution gm of oxygen equiv $= \frac{8aV}{10}$ gm of oxygen

20 ml of water contains equiv $\frac{8aV}{10}$ mg of Oxygen

1000 ml of water contains of equiv $= \frac{8aV}{10} \times \frac{1000}{20}$ mg/L of Oxygen

= "A" mg/L of Oxygen

= "A" ppm of Oxygen

Conclusion:

The amount of dissolved oxygen present in the supplied sample of water is..... ppm.

Experiment :6

Aim of the Experiment : Determination of Percentage of Available Chlorine in a supplied bleaching powder solution.

Apparatus Required:

Burette, Pipette, Conical Flask, Watch Glass, Measuring Cylinder, Wash Bottle

Chemicals Required:

Bleaching Powder Solution, KI (10%) Solution, Dilute Acetic Acid, (N/10) Sodium thiosulphate Solution, Starch Indicator

Theory :

Bleaching powder is a good disinfectant which has added to water for killing disease-causing bacteria. The germicide action of bleaching powder is due to hypochlorous acid and nascent oxygen liberated. It is powerful oxidizing agent and hygroscopic in nature.



The quality of bleaching powder is based on the quantity of available chlorine presenting in it. Percentage of available chlorine may be defined as the amount of chlorine available for oxidation per hundred parts of bleaching powder. So, more is percentage of available chlorine, better is the quality of the bleaching powder. Although, bleaching powder contains about 56% of chlorine but a good quality of bleaching powder contains only 35-38% of available chlorine. Because on standing bleaching powder undergoes slow decomposition and gets converted into calcium chloride and calcium chlorate. Hence percentage of available chlorine decreases on passage of time.

Principle:

The chlorine, which is liberated in water made to react with KI so that iodine is liberated. This liberated iodine is titrated against a standard solution of sodium thiosulphate using starch as indicator.

Chemical Reactions:

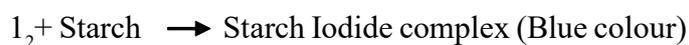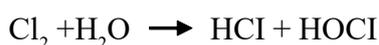Procedure:

Pipette out 10 ml of supplied bleaching powder solution into a 250 ml of conical flask. Add 15 ml of 10% KI solution followed by 10 ml of dilute acetic acid solution. Cover the conical flask with watch glass and set for 5 minutes in dark. After 5 minutes, titrate the liberated iodine with standard hypo-solution to get a faint yellow colour. Add 2-3 drops of starch indicator and continue the process of titration to get a colorless solution, which is the end point. Note the titrate value. The above process is repeated for another 3 times to get the concordant readings.

Tabulations :

No of Observations	Volume Supplied CaCOCl ₂ solution of	Initial Burette Reading in ml	Final Burette Reading in ml	Difference (ml)	Remark
1					Rough
2					Concordant
3					
4					

Calculations :

Concentration of Bleaching powder solution = 10gm / L

Titrate value = Xml

We know that :

1000 ml of 1 N Na₂S₂O₃ = 35.5 gm of chlorine

1 ml of 1N Na₂S₂O₃ = 0.035 gm of chlorine

1ml of “a” N/10 Na₂S₂O₃ ≈ 0.035a gm of chlorine

X ml of “a” N/10 Na₂S₂O₃ ≈ 0.035aX gm of chlorine

Now 10 ml of supplied bleaching powder solution contains 0.035aX gm of chlorine

1000 ml of supplied bleaching powder solution contains $\frac{0.035aX}{25} \times 1000 = A$ gm of chlorine

% of Available Chlorine :

10 gm of CaOCl₂ solution contains A gm of available chlorine

100 gm of CaOCl₂ solution contains 10A gm of available chlorine

Conclusion :

Supplied bleaching powder solution contains of available chlorine.

Experiment :7

Aim of the Experiment : Estimation of amount of calcium in lime stone.

Apparatus Required : Burette, Pipette, Conical Flask

Chemicals Required:

KMnO₄ Solution, Standard oxalic acid (N/10), 8% Ammonium oxalate solution, NH₄Cl, Ammonia

Theory : Lime stone essentially consists of calcium carbonate but it is generally associated with small quantities of magnesium carbonate. Dolomite is an equimolecular compound of calcium carbonate and magnesium carbonate. These ores are generally contaminated with minor amounts of oxides of iron, aluminium, silicon (free or combined) and organic matters. The minerals have large number of industrial applications and the quality of the mineral as well as its suitability for a particular commercial application depends upon the amounts of CaCO₃, MgCO₃ and the other constituents present. Hence, complete analysis is often necessary.



In this method, the lime stone powder is dissolved in hydrochloric acid and the calcium present in the solution is precipitated as calcium oxalate using oxalic acid or ammonium oxalate in presence of ammonia. The precipitated calcium oxalate after washing is treated with dilute sulfuric acid and the oxalic acid liberated is titrated with a standard KMnO₄ solution. From the volume of KMnO₄ required for the titration, the amount of Ca in the ore can be calculated:

$$1 \text{ ml of } 1 \text{ N KMnO}_4 = 0.020 \text{ gm of Ca}$$

$$1 \text{ ml of } 1 \text{ N KMnO}_4 = 0.028 \text{ gm CaO}$$

$$1 \text{ ml of } 1 \text{ N KMnO}_4 = 0.050 \text{ gm of CaCO}_3$$

Chemical Reactions :

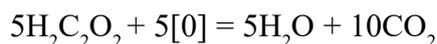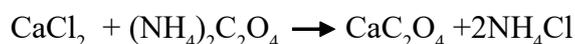Procedure :

Weight accurately 1 gm of the lime stone in to a beaker and treat it with a little dilute HCl to dissolve the ore. Take care to see that there is no loss of the solution due to quick evolution of CO₂ by covering the beaker with a watch glass. After all the ore is dissolved, wash the watch glass with water in to the same beaker. Render the solution alkaline by adding ammonia. Add 1 gm of NH₄Cl and stir the solution to

dissolve it. Continue the stirring and add excess of 8% ammonium oxalate solution. Boil the solution and allow the calcium oxalate precipitate to settle down. Decant the supernatant liquid through a filter paper. Wash the precipitate several times with hot distilled water containing a little of ammonium hydroxide until the filtrate is free from chloride and oxalate ions (by testing with AgNO_3 for chloride and mild hot and dilute acidified KMnO_4 for oxalate). Preserve this filtrate for the determination of magnesium. if necessary, as magnesium pyrophosphate.

Place the beaker containing the precipitate under the funnel and add 20 to 25 ml of 5N H_2SO_4 in to the filter paper to dissolve the calcium oxalate present in the filter paper. The acid also dissolves the calcium oxalate precipitate in the beaker. Wash the filter paper with distilled water and collect these washing also in to the same beaker. Transfer the solution quantitatively in to a 250 ml measuring flask, make it up to the mark and shake the solution thoroughly. Take 25 ml of the solution in to a conical flask, add 10 ml of dil. H_2SO_4 , heat the solution to about 70°C and titrate against the standard KMnO_4 solution (which in turn was standardised with standard oxalic acid solution first).

No of Observations	Volume of Oxalic Solution	I.B.R. (ml)	F.B.R (ml)	Difference in ml	Remarks
1					Rough
2					Concordant
3					
4					

Calculations :

Weight of ore taken = 1 gm

Volume of the oxalic acid solution made up = 250 ml

Volume of oxalic acid taken for each titration = 25 ml

Volume of KMnO_4 solution rundown = x ml

Standardisation of KMnO_4 solution :

25 ml of N/10 oxalic acid = 25 ml of N/10 KMnO_4 solution

Normality of KMnO_4 = N/10

Now, we know that

1 ml of 1 N KMnO_4 = 0.02 gm of Ca

19.3 ml N/10 KMnO_4 = $0.02 \times X \times 1/10$ gm of Ca

Hence, 25 ml of the solution contains $0.02 \times x$ gm of Ca

250 ml of the solution contains $0.02 \times x$ gm of Ca

Now, 1 gm of the ore contains $0.02 \times x$ gm of Ca

100 gm of the ore contains $0.02 \times X \times 100$ gm

= x gm of Ca

% of Ca present in the ore = $38.6\% = 96.5\% \text{CaCO}_3$

Conclusion:

The supplied lime stone contains $x\%$ of Ca.

Experiment :8

Aim of the Experiment : Determination of flash point of an oil using Pensky-Martens flash point apparatus

Apparatus Required:

Pensky-Marten's Closed Cup Flash Point Apparatus, Thermometer, Measuring Cylinder, LV Beaker

Chemical Required:

Supplied Oil Sample

Theory:

The flash point of an oil is defined as the minimum temperature at which it gives vapours, which when mixed with air form an inflammable mixture and gives a flash of light on the application of a pilot flame. This gives information regarding the volatility of oil necessary to prevent fire hazards.

The fire point is the lowest temperature at which sufficient vapours are grown up by an oil, so that it burns continuously for at least 5 seconds on being lighted by a pilot flame.

The flash point gives a rough indication of the tendency of the oil to volatilise and therefore is useful in process control for maintaining the uniformity of the product. The fire point of an oil is about to 40° higher than its flash point. A lubricant should have its flash point at least above the temperature at which it is to be used. This provides protection against fire hazards during use, transport and storage of lubricants. Both flash and fire points can be determined by using Pensky-Martens closed cup apparatus.

Procedure :

Take the supplied oil sample in the cup provided in the apparatus. Close the cup with the cover and place the thermometer in its position. Adjust the electrical heating system so that the rise in temperature is about 1 °C in 7 minutes. Keep the pilot flame burning and continue hosting with constant stirring. Take the reading at every 5 °C interval to obtain the approximate flash point. Allow the oil to cool with stirring by 5 °C and take the readings with every 1 °C rise in temperature till the flash point is obtained. Once the flash occurs, allow the temperature to rise by about 5 °C more and discontinue the heating. Allow the oil to cool with continuous stirring and take the readings with every 1 °C fall of temperature till no flash occurs. The minimum temperature at which the flash occurs gives the flash point of the oil.

Tabulation: 01 (Observation for approximate flash point)

No. of Observations	Temperature increase by (5 °C)	Remarks
1		
2		
3		
4		
5		
6		
7		
8		
9		
10		

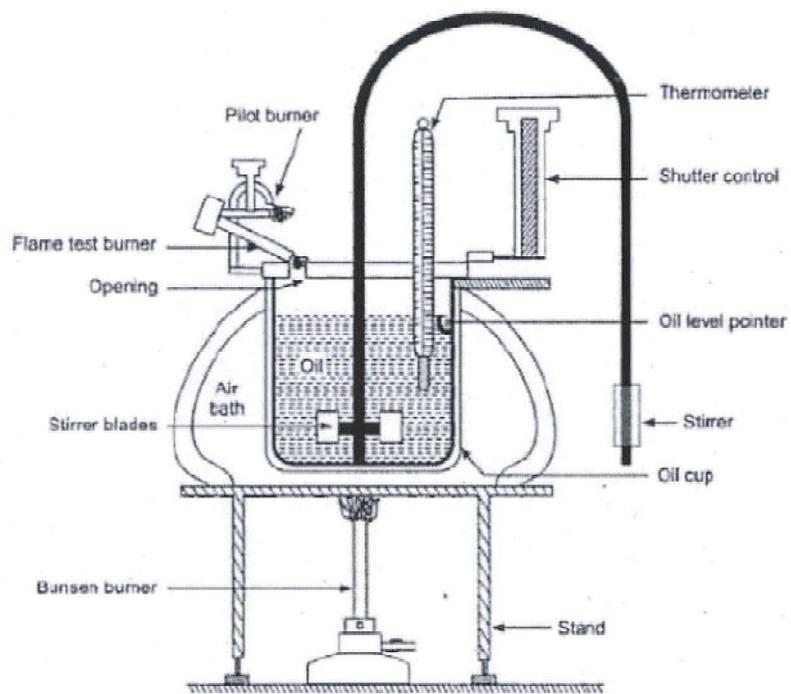
Figure 01: Pensky Marten's Flash Point Apparatus

Tabulation: 02 (Determination of flash point)

No. of Observations	Temperature Increase by (5 °C)	Remarks	No. of Observations	Temperature Increase by (5 °C)	Remarks
1			1		
2			2		
3			3		
4			4		
5			5		
6			6		
7			7		
8			8		
9			9		
10			10		

Conclusion :

The flash point of the given oil is _____ °C = _____ °F.

The fire point of the given oil is _____ °C = _____ °F.

Experiment : 9

Aim of the Experiment : Determination of viscosity of a lubricating oil by using Redwood viscometer.

Apparatus Required :

Redwood Viscometer, Thermometer, Stopwatch

Chemical Required :

Given Lubricating Oil

Theory :

The **Viscosity** is the property of a fluid by virtue of which it offers resistance to its flow. It is thus a measure of the internal resistance in the fluid and its determination is usually done in case of the lubricating oils. It is the single most important property of any lubricating oil, since it is the main factor that determines the operating characteristics of a lubricant.

The **Absolute Viscosity** of an oil at any given temperature is determined by measuring the rate of flow of the oil through a capillary tube. The **Kinematic Viscosity** of lubricating oils is generally determined by measuring the time taken for a given quantity of the oil to flow under its own weight through a capillary of standard dimension under standard conditions. The instrument used for the measurement of viscosity of the oil is called viscometer. **Redwood Viscometer** is commonly used for measurement of viscosity of lubricating oils, in empirical unit and not in absolute unit.

Absolute and Kinematic Viscosities of a fluid can be determined by **Redwood viscometer**, according to the formula

$$\text{Absolute Viscosity} = \left(\alpha t - \frac{\beta}{t} \right) \rho$$

$$\text{Kinematic Viscosity (V)} = \frac{\text{Absolute Viscosity}}{\rho} = \alpha t - \frac{\beta}{t}$$

Where,

V = Kinematic viscosity of the oil in centistokes.

T = Time of flow in seconds for a fixed quantity of the oil flowing through capillary

ρ = Density of the oil

α and β are instrument constants.

The value of α and β are in given below:

Time of flow in seconds for 50 ml	α	β
34 - 100	0.00260	0.00260
Over 100	0.00247	0.00247

The values are calculated in centistokes unit.

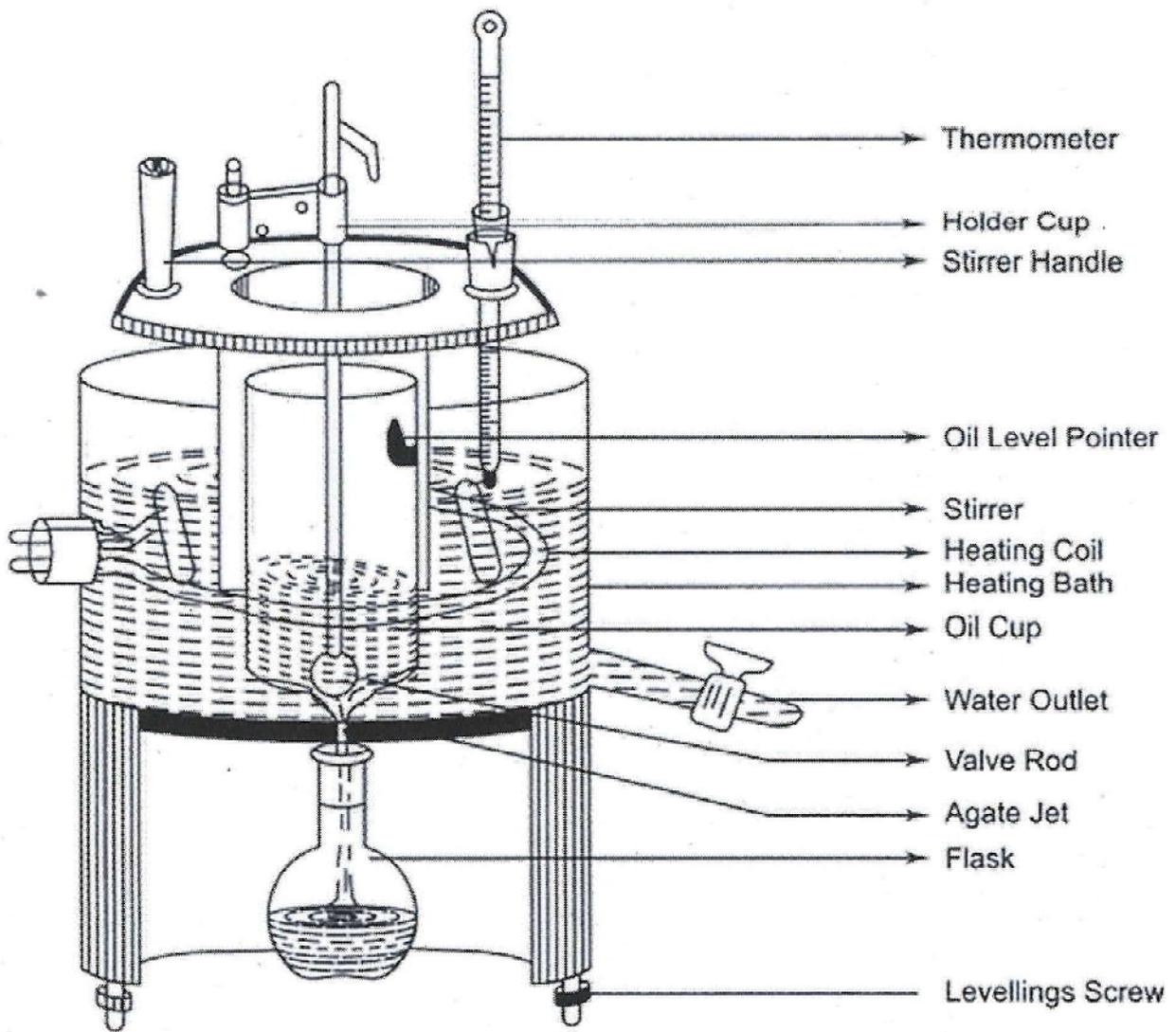
Figure 01: *Redwood Viscometer.*

Procedure :

Before taking the oil in the oil cup ensure that the oil cup is clean and dry and the capillary has no obstructions for the flow of the oil. Close the capillary and fill the cup with oil sample up to the level indicator. The heating jacket is filled with water, which can be heated to rise the temperature of the oil to different values at which measurements are to be made. Make the first measurement when the oil is at room temperature by allowing the oil to flow through the capillary and record the time in seconds when exactly 50 ml of the oil is collected in the measuring Kohlrusch flask placed at the bottom. Close the capillary, refill the oil cup by the same oil up to the levelling mark, increase the temperature of the oil by about 10-15 °C (not less than 10 °C) by heating and record the temperature and the corresponding time in seconds. when exactly 50 ml of the oil is collected in the measuring Kohlrusch flask. Repeat the measurements at five different temperatures. Calculate the kinematic viscosity of the oil at different temperatures and plot a graph between temperature and kinematic viscosity of the oil. Report the mean value as Redwood viscosity in seconds.

Tabulation :

No of Observations	Initial temperature (°C)	Final temperature (°C)	Mean temperature (°C)	Time (Second)	Value of α	Value of β	Kinematic Viscosity (centistokes)

Graph : Plot the graph Temperature versus Kinematic Viscosity

Conclusion :

The downward sloping curve obtain by plotting a graph of temperature versus kinematic viscosity proves that the above are inversely related to each other i.e. with increase in temperature, viscosity decreases and vice-versa and the viscosity of the supplied oil is found to be Redwood seconds.

Experiment : 10


Aim of the experiment : Preparation of Aspirin (acetyl salicylic acid)

Apparatus required : Conical flask, Beaker (100 mL), Buchner funnel

Chemicals required : Salicylic acid, Acetic anhydride, Conc. H_2SO_4

Theory:

Salicylic acid upon acetylation yields acetyl salicylic acid or aspirin. Acetylation proceeds rapidly with acetic anhydride in presence of a little conc. H_2SO_4 as catalyst.

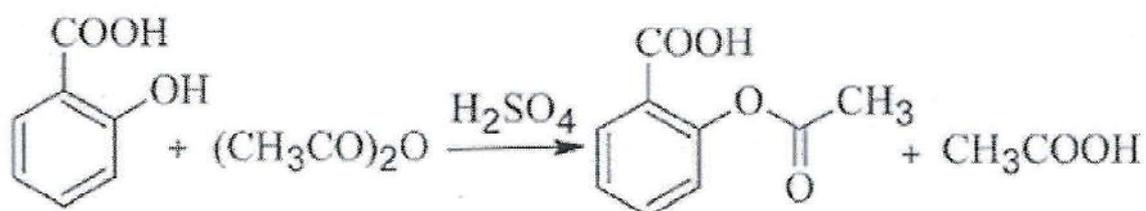Procedure :

Take 5 g of salicylic acid and 7 mL of acetic anhydride in a 100 mL. conical flask. To it add 2-3 drops of conc. H_2SO_4 ; ; shake and then heat the flask on a water bath keeping temperature between 60-70 °C for 20 minutes. Cool the mixture in an ice bath with stirring; add about 20 mL of the ice cold water to decompose the excess of acetic anhydride. Pour the contents of the conical flask into a 250 mL beaker containing 50 g of crushed ice, stir until white solid appears. Scratch the side wall of the flask with a glass rod or stainless steel spatula, if white solid does not appear. Filter through a Buchner funnel, dry by pressing into the folds of filter paper. Recrystallize the crude acetyl salicylic acid from a mixture of equal volumes of acetic acid and water.

Conclusion:

Yield-5g, M.P=136-137 °C

Acetyl acetic acid decomposes when heated and does not possess a true clearly defined melting point

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